
scala®

Digitales Fieberthermometer SC 35 T, SC 37 T, KD-132

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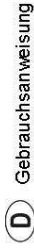CE 0197

SCALA Electronic GmbH
Ruhlfelder Str. 95
14532 Stahnsdorf

Digitales Fieberthermometer SC35 T, SC 37 T, KD-132

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Clinical thermometer SC35T, SC 37 T, KD-132

Special features:

The SC35T, SC 37 T, KD-132, is a maximum thermometer, i.e. the highest temperature is displayed after only 60 seconds. The most recent measured value is always saved.

Safety information:

- The SC 37 T is suitable only for measuring human body temperatures.
- The thermometer contains small parts (batteries, etc) which could be swallowed by small children. Please keep the thermometer and its batteries out of the reach of small children. If the battery or other parts are swallowed you should seek medical help immediately.
- The enclosed batteries may not be recharged or reactivated by other methods, may not be dismantled, thrown into fires or short-circuited. Batteries must be disposed of properly. Suitable methods for the disposal of batteries are given in the instructions for use. Please refer to the disposal instructions for batteries in the place of purchase or disposal in communal collection points.
- Protect the SC 37 T from direct sunshine and temperatures above 50°C.
- Never drop the SC 37 T, do not bend or open it (with the exception of the battery compartment for battery replacement).

Operating manual:

Press the I/O button to switch on. A short signal tone shows that the device is in the "switched on" condition. A segment test takes place simultaneously (checks all display segments and carries out an internal function test (see under "Function test").

After this, if the temperature measured by the temperature sensor is less than 32°C, the display shows an "L" and "°C". After one and one second, the "°C" starts to blink. The "L" stands for low (less than 32°C).

The current temperature is continuously displayed during measurement, and the "°C" blinks during the entire measurement time. The measurement is completed after around 60 seconds.

If the measurement is taken orally or rectally, the approximate final temperature is shown after about 60 seconds and changes only marginally after this time. A measurement period of between five and 10 minutes must be observed for the underarm measurement method. The signal tone can be sounded several times during this period. Since the contact between the temperature sensor and the armpit is not as good

as those for oral and rectal measurements, the temperature rise is slower which causes this effect.

After the measurement is completed, the temperature value measured remains in the display until the SC 37 T is switched off. To turn the device off after use, press the I/O button briefly. The signal tone will again be sounded. If the device is not switched off, a signal tone is sounded after around 8 minutes and the device then switches off automatically.

Display of last measured value:

If you keep the I/O button pressed for longer than three seconds, the last measured value will be shown in the display. As soon as you release the I/O button, the device is ready to take measurements.

Function test:

The thermometer is checked for correct functioning automatically every time it is switched on. All segments are displayed, and then L.C. If a failure (measurement inaccuracy) is recognised, ERR (Error) is shown in the display.

Measured value deviation:

The measured temperatures will always deviate, depending on the measurement method, from the actual temperature. The deviation in temperature distribution on the human body. You could, for example, measure a rectal temperature of 36.6°C, an oral temperature of 36.1°C and an underarm value of 35.6°C.

Measurement locations:

In the mouth (Oral):
Small temperature variations in the mouth between measurements are normal. These are caused by varying temperature distribution within the mouth. Try to place the thermometer in one of the two heat pockets under the tongue, to the left and right of the tongue root, and ensure the tongue root must have good contact with the tissue. The subject should close the mouth and breathe normally through the nose so that the measured result is not affected by breathed air.
In the anus (Rectal):
This is the safest method. It is particularly suitable for infants and small children. Guide the thermometer point carefully 2 - 3 cm into the anus.
Under the arm (Axillary):

From a medical viewpoint this type of measurement leads to inaccurate values. It is therefore unsuitable for carrying out exact measurements.

Measuring period: 5 to 10 minutes.

Measuring the base/body temperature:
The device can also be comfortably and practically used for the measurement of fertile or unfertile days in the female cycle using the temperature measurement method.

Changing the battery:

As soon as the symbol appears in the display, the battery has run out and requires replacement. Always replace the discharged battery with one of the same type. To do this, proceed as follows:
• Pull the cover off the battery compartment. The battery can now be accessed.
• Put the new battery in the battery compartment with the correct polarity (SR41 in. Please make sure it is the right way round. You can purchase batteries from jewellers or electrical shops.

Measurement checks:

Checking the accuracy of the thermometer is not necessary in most cases. The German Medical Product User regulations of 29th June 1998 stipulate a measurement check every two years for devices used in medical practices. If you require a measurement check (at your expense), please send the thermometer to our customer service address. Please make sure it is well-packaged and that the correct postage stamps are applied.

Cleaning and disinfection:

The thermometer is water-tight, and can therefore be fully immersed in one of the cleaning fluids listed below for cleaning.

Name of fluid	Manufacturer	Application
Edisonite	Merz GmbH & Co.	in accordance with manufacturer's instructions
Edisonite (Super)	Merz GmbH & Co.	in accordance with manufacturer's instructions
Gigaset PF	Schülke & Mayr GmbH	in accordance with manufacturer's instructions

We recommend the use of Mucocit®F2000, or Helipus®Plus, to disinfect the device during use in medical facilities. Both are effective against all hospital-typical pathogens and can be purchased from any dispensing chemist.

Guarantee:

We give a guarantee period of 2 years from the date of purchase. We guarantee perfect design of this article. We will rectify all material or manufacturing errors during the guarantee period. If, contrary to expectations, defects are found, please send the device back to our customer service department. Please enclose the

proof of purchase and a description of the defect so that the problem can be rectified as soon as possible. The guarantee excludes damages caused by improper handling, related to the battery or to the packaging.

Referenced norms:

The device has been certified in accordance with DIRECTIVE 93/42/EU of 14th June 1993 relating to medical products.
The clinical thermometer SC35 T, SC 37 T, KD-132 complies with the requirements of European norm EN 12470-3:2000.

Disposal note

The disposal of batteries is regulated. Spent batteries are collected at collection points provided by battery dealers and the communities. Dispose of the thermometer monitor as prescribed, i.e., at a collection point for electrical scrap.

Specifications:

- maximum thermometer
- 32.0°C to 43.9°C
- Temperature below 32.0°C: display "L" for low.
- Temperature above 43.9°C: display "H" for high.
- Measurement accuracy: 0.1°C between 35.5°C and 42.0°C, at a surrounding temperature of 18°C to 28°C.
- Automatic internal check of test value 37°C.
- Increase or deviation >0.1°C, "Err" will be displayed.
- 3-dot liquid-crystal (LCD display).
- 1 no. LR41 or SR41.
- The "D" shows when the battery is discharged.
- Storage temperature: 0°C to 50°C

Explanation of symbols

Symbol	Manufacturer	Product type BF	Disposal notes
	SCALA Electronic GmbH	Medical device	Read the instruction manual
	SCALA Electronic GmbH	Medical device	Read the instruction manual
	SCALA Electronic GmbH	Medical device	Read the instruction manual
	SCALA Electronic GmbH	Medical device	Read the instruction manual
	SCALA Electronic GmbH	Medical device	Read the instruction manual

Service client:
SCALA Electronic GmbH
Ruhlsdorfer Str. 95
14532 Stahnsdorf
Germany

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Senso elektronik, Stahnsdorf GmbH